


Battery Operated Handpiece

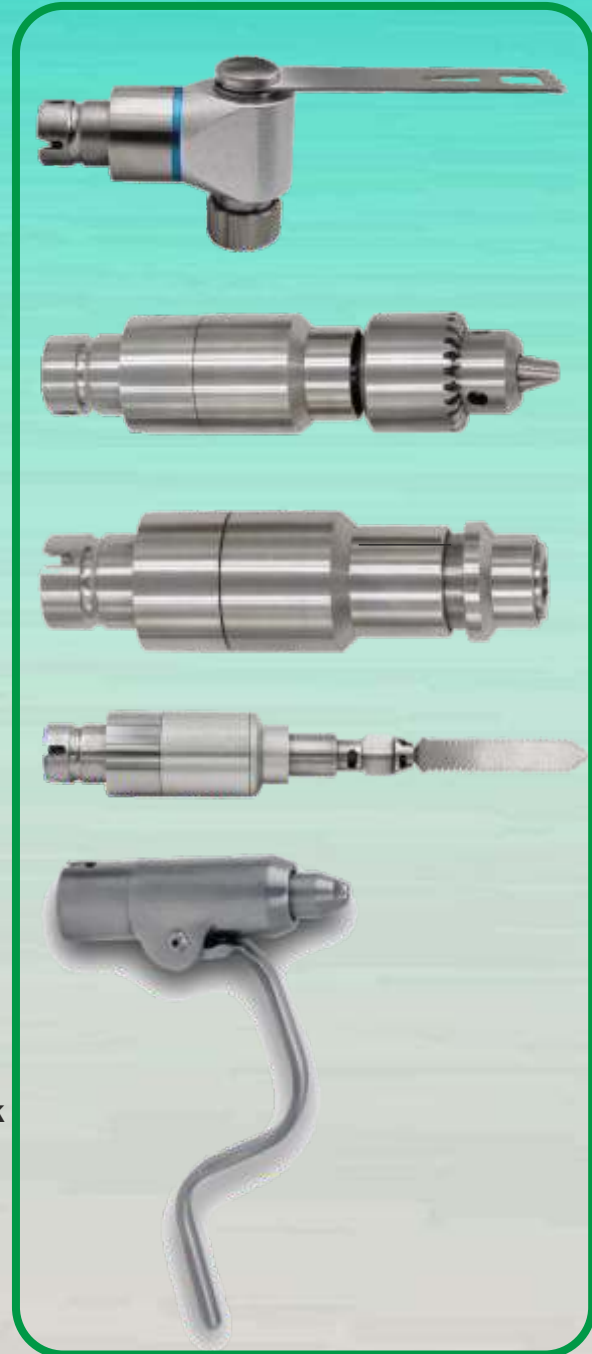

FROM WHEELS TO HEELS

PRODUCT APPLICATION

1. Apply to all kinds of brands of intramedullary nail, k-wire, trauma and joint operations.
2. The handpiece has perfect shape of ergonomic design, users feel comfortable to handle.
3. The handpiece is integrated design, and it has good sealing. The whole unit can be autoclavable.
4. The rotating lock is securer and more durable.
5. The handpiece uses special material, and it has good wear resistance and water resistance. Better than medical stainless steel.
6. High-quality motors. Low loss, low noise, low wear and longer service life.
7. Long-lasting and consistent high-capacity batteries.
8. Adaptors: Sagittal Saw, Reciprocating Saw, Drill Chuck, Reamer Chuck, K wire Driver (0.8mm-1.8mm), K Wire Driver (2.0mm-3.2mm)

2. Specification

| Specification | | Standard Configuration | |
|-------------------------|--------------------|----------------------------|-------------------|
| Input Voltage | 110V-240V/50Hz | Handpiece | 1 |
| Diameter of drill chuck | 0-8 mm | Adaptors | Can be customized |
| Output Power | ≥ 120W | Battery | 2 |
| Battery Capacity | 1800mA/h | Key | 1 |
| Rated Speed | 0-1200rpm/20000rmp | Charger | 1 |
| Rating torque | 6.3Nm | Aseptic Battery Tranfering | 1 |
| Charging Voltage | 14.4V | Operating Manual | 1 |
| Sterilizing Temperature | 135-137°C | Guarantee Card | 1 |
| Noise | ≤60db | Alluminum Case | 1 |

3. Warranty

The Handpiece guarantees for 1 Year;

The battery and other fast consumption items guarantee for 3 months.

4. Notes

- *Sterilize temperature is 135-137°C.*
- *After the operation, open the cover of the battery compartment and take down the battery pack.*
- *The battery of this product must not be disinfected under the conditions of high temperature and high pressure, not be kept in the machine for a long time.*
- *No water in battery space after sterilization.*
- *If you find the abnormal sound from the machine or the overheating of the machine in use, please stop using, and contact the manufacturer or distributor immediately, and then send the product to the manufacturer for repair.*




↑
ACE-701

Multi-functional handpiece




ACE-702

Sagittal Saw Attachment

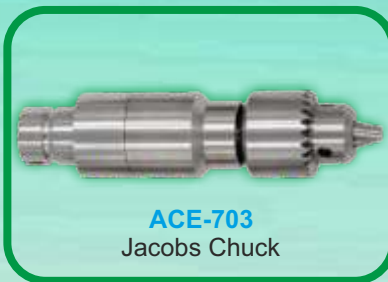

ACE-703

Jacobs Chuck

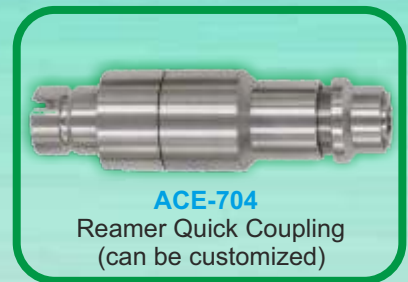

ACE-704

Reamer Quick Coupling
(can be customized)



ACE-705

Reciprocating Saw Attachment



ACE-706

Pin & Wire Driver
0.8-1.8 mm
(For cannulated
handpiece)




ACE-707

Pin & Wire Driver
2.0-3.2 mm
(For cannulated
handpiece)



ACE-708

Cannulated Drill Attachment

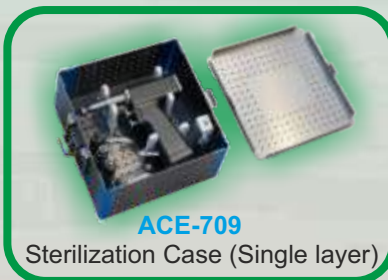

ACE-709

Sterilization Case (Single layer)



ACE-710

Aseptic Battery
Transfer
Kit+C21:C22



ACE-711

Key



ACE-712



NiMH Battery



ACE-713

Charger

www.orthocare.in
www.orthoimplants.com

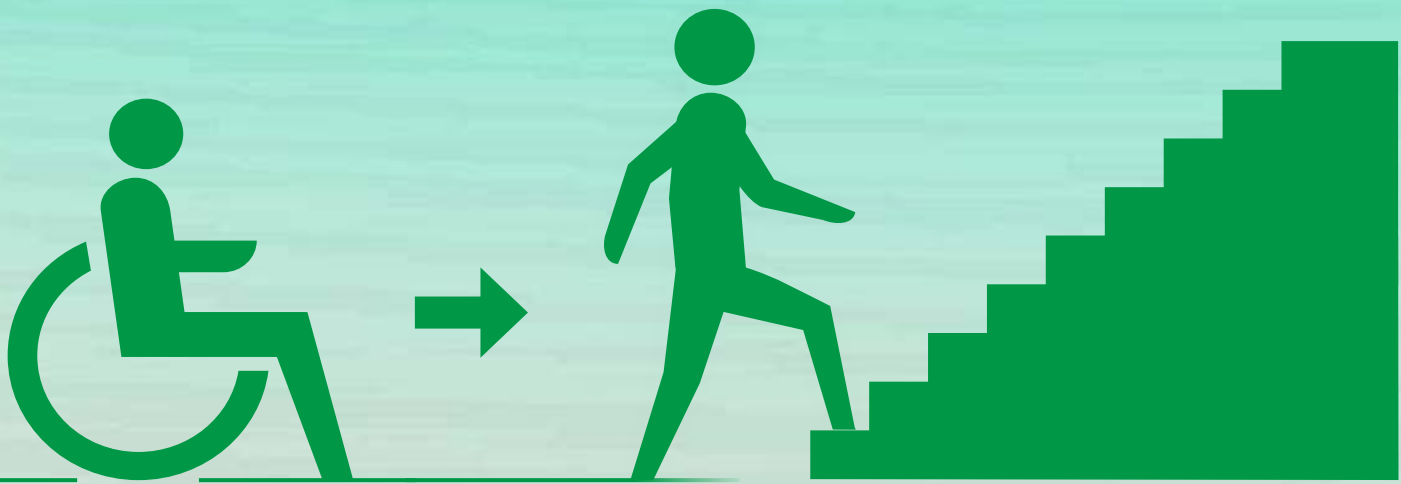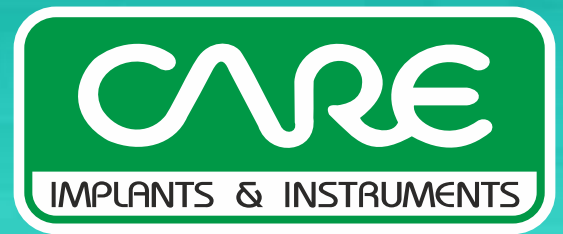

FROM **WHEELS** TO **HEELS**

ORTHO CARE

H.O.: 3/10 Kirti Nagar Industrial Area, New Delhi – 110015, INDIA
Phone : +91-11-41424511 | 41424522 | 45092044 Fax : +91-11-47510708
E-mail : orthocareindia2009@gmail.com Web : www.orthocare.in | www.orthoimplants.com
Factory : A-31, Sector A-4, Tronica City Loni, Ghaziabad (U.P.) INDIA

USFDA 510K CLEARED PRODUCTS
CE 1023 CERTIFIED PRODUCTS
ISO 13485:2012 CERTIFICATION
INDIAN FDA LICENSE NO. 05/SCP/2011